

iTeraCare™



iTera-*Pro*

Prife International Sdn Bhd

A-GR-2, Menara Uncang Emas, No 85 Jalan Loke Yew, Cheras 55200 Kuala Lumpur, Malaysia



www.i-teracare.com



info@i-teracare.com


iTera-Pro

Dear User :

Thank you for purchasing iTera-Pro device.
Please read the Manual carefully before using.

Product Description :

- **Product appearance** (see diagram)
- **Product model** : DS-808
- **Voltage** : 110V
- **Frequency** : 50Hz
- **Power** : 1200W
- **Sound** : $\leq 65\text{dB}$

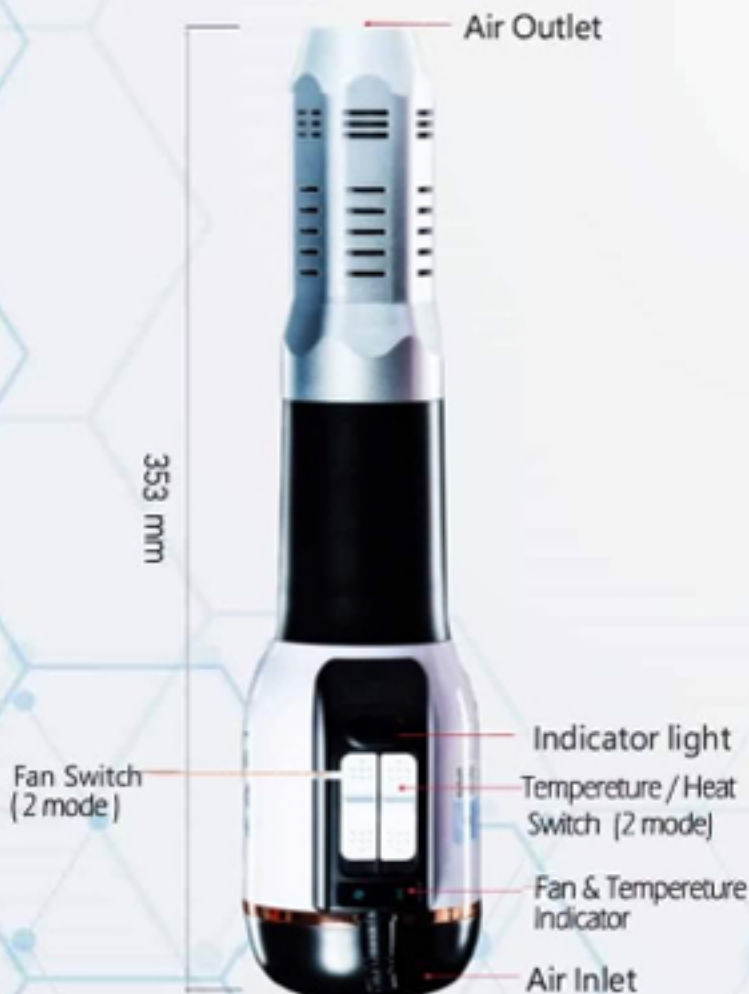

Warranty description :

One year warranty (from date of purchase)

1. If the product is found defective as a result of manufacturing defects in material and workmanship during the warrant period, please contact our customer service.
2. Company will provide repair and maintenance service after warranty period where cost will be incurred to the owner.

Warranty will be void under the following conditions :

1. Exceeded warranty period.
2. Physical damage due to not following user manual provided.
3. Other circumstances that void the warranty :
 - ▶ Device mishandling.
 - ▶ Defects or damage caused by unauthorized repair.
 - ▶ Defects or damage due to fall and impact.

Environmental Protection device :

- This device has Environmental protection usage for 10 years
- The Environmentally Protection assures that the device is not harmful to humans and environmentally friendly.

CAUTIONS :

1. Read all instructions carefully before using.
2. Handle with care. The device is fitted with a crystal tube, avoid knocking or dropping the device.
3. Do not block the air inlet and outlet to prevent overheating.
4. Remove all metal accessories before using. Remove contact lenses before using on the face.
5. To protect against electrical shock, do not immerse cord, plugs, device body in water or other liquids.
6. Unplug from mains when the device is not in use.
7. Do not use this device at outdoors. This device is intended for household use only. Do not use the device for other than its intended use.
8. Do not use it for too long period of time. This device is fitted with a resettable thermal link to protect the motor from overheating. Should the motor overheat, the thermal link will disconnect the power supply and will automatically re-connect the power supply after the motor cools down.
9. If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
10. This device is not intended for use by persons (including children) with reduce physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety. Children should be supervised to ensure that they do not play with the device.
11. Never leave the device unattended when it is switched on.
12. Do not switch on the device if it appears to be faulty in anyway. Please contact and send to nearby service center.